

Simplifies the administration of medical oxygen to patients

Your patient is scheduled for a CAT scan. You have a cylinder available but can't find a regulator or a wrench. The patient's test is delayed, and radiology has to reschedule. Your facility suffers lost time, money, and quality of patient care.

Now you can prevent oxygen-related transport delays – with the Grab 'n Go™ portable medical oxygen system – a unique technology that gives you USP oxygen at your fingertips. The regulator and contents gauge are built right into the cylinder.



No-risk trial

The best way to determine whether the *Grab 'n Go* oxygen system is right for your facility is to try it absolutely free! For a limited time, you can use up to 20 *Grab 'n Go* oxygen units in your facility, for up to 3 months, at no cost.

For more information, call 800-299-7977 ext. 6961, contact your Praxair representative, or visit www.praxair.com/healthcare

Why choose the Grab 'n Go™ oxygen system?


- Reduced equipment costs The system does not require a separate regulator to be attached. This eliminates the need to buy, inventory, or repair regulators for oxygen cylinders.
- Less labor Just pick up a Grab 'n Go oxygen system and turn the dial to the marked setting. You don't have to find a wrench and attach a regulator to the oxygen cylinder prior to use.
- Simplified training One-knob operation makes the Grab 'n Go system the simplest medical oxygen equipment in the industry today. And we offer a free training video and CD to further reduce the learning curve.
- Increased safety High impact plastic cover protects the regulator in the event of a drop or fall.
- Eliminates transport delays You don't have to spend time searching for a regulator to take your patient to in-house diagnostic procedures, or calling the respiratory therapy department to straighten things out. This prevents delays and increases diagnostic throughput.
- Meets emergency requirements An inventory of Grab 'n Go oxygen systems gives you the supply of patientready oxygen needed to evacuate your facility in case of emergency.

An all-in-one system

Praxair offers the complete Grab 'n Go™ oxygen system as a single integrated unit. There are no parts to assemble or keep in inventory, no tools required for use.

Each Grab 'n Go portable medical oxygen system includes:

- Standard E cylinder of medical-grade USP oxygen special lightweight construction results in a total system weight of just 12 pounds.
- Attached solid brass regulator constructed from brass, the preferred material for oxygen regulator applications, and eliminates the need to attach regulators or change-out cylinders. When a Grab 'n Go oxygen unit is empty, just return the entire system for a replacement.
- Built-in contents gauge continuously indicates pressure to show whether the cylinder is full or empty,
- Convenient carrying handle ergonomically designed for easy handling and transport.
- High-impact plastic housing protects the regulator from shocks and damage.
- Snap-set flow rate control permits flow rates to be easily set in the half-liter increments frequently prescribed by doctors.

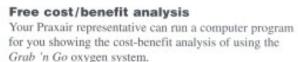

The regulator and maintenance savings alone can often offset the entire cost of the *Grab 'n Go* program. There are additional savings including increased diagnostic revenues achieved through more timely patient transport, a shorter learning curve for your staff, and improved safety.

Why hospitals and nursing homes prefer Praxair

Thousands of *Grab 'n Go* portable medical oxygen systems are in use at hundreds of health care facilities throughout the U.S. and Canada. These portable units are available from more than 300 Praxair locations in North America.



PRAXAIR

© 2001, Praxair Technology, Inc. All rights reserved.

Praxair, Inc. 39 Old Ridgebury Road Danbury, CT 06810-5113 USA

www.praxair.com/healthcare

Tel: 1-800-299-7977 ext. 6961

Praxair, the Flowing Airstream design, Grab 'n Go and Grab 'n Go III are either trademarks or registered trademarks of Praxair Technology, Inc. in the United States and/or other countries. The information contained herein is offered for use by technically qualified personnel at their discretion and risk without warranty of any kind.

Manufactured for Praxair by Western Medica, Westlake, OH 44145

P-8578 01/01 12M